

800 Series Benchtop pH/Conductivity Meter User Manual

- ☐ PH800 Lab pH Meter
- ☐ EC800 Lab Conductivity Meter
- ☐ PC800 Lab pH/Conductivity Meter



IP54

Table of Contents

1. Introduction	4
2. Configuration	6
3. Technical Specifications	7
4. Instrument Description	8
5. pH Measurement.....	14
6. mV Measurement.....	24
7. Conductivity Measurement	25
8. Parameter Setting	35
9. USB Data Communication.....	42
10. Recommended pH Electrodes for Different Applications.....	46
11. Recommended Conductivity Electrodes for Different Applications ..	47
12. Replacement Parts and Accessories.....	48
13. Warranty	49
14. Appendix I : Table of Parameter setting and factory default setting ..	50
15. Appendix II: Icons and Abbreviation	51
16. Appendix III : Table of Self-diagnosis Symbol	53

Special Notes

- Do NOT pull out the power plug until the meter is turned off.
- Don't pull out USB plug or press  button to turn off the meter directly when the meter is connecting with PC. Please press "Exit" button first (in the PC interface) to quit PC Link Software, then press  button to turn off the meter, and pull out USB plug after that.
- There is a CR2032 3V lithium battery in the meter's circuit board, see the picture on the right. It's power source for meter's inner clock. When the voltage is lower than 2.6V after long time use, a blue screen with on display values and icons may appear when turn on the meter. At that moment, please open the meter and replace a new CR2032 3V lithium battery.



⚠ CAUTION



Chemical exposure hazard. Obey laboratory safety procedures and wear all of the personal protective equipment appropriate to the chemicals that are handled. Refer to the current safety data sheets (MSDS/SDS) for safety protocols by scanning the QR code on the right.



1. Introduction

Thank you for purchasing 800 Series Benchtop pH/Conductivity Meters. Before using the product, please read this manual carefully to help you properly use and maintain the product.

1.1 Measuring Parameters

Measuring Parameters	PH800	EC800	PC800
pH/mV	√		√
Conductivity/TDS/Salinity/Resistivity		√	√
Temperature	√	√	√

Note: PC800 pH/Conductivity Meter does NOT support simultaneous measurement of pH and conductivity.

1.2 Features and Functions

- The built-in microprocessor chip enables advanced functions such as auto calibration, auto temperature compensation, auto electrode recognition, parameter setting, self-diagnosis, calibration reminder, calibration time check, auto power-off, low-battery reminder, etc.

- GLP data management, real-time clock display. Manual or Auto timing data storage. USB data communication.
- The meter adopts advanced digital processing technology, intelligently improves the response time and accuracy of the measurements. Stable reading and auto lock display mode are available for choice.
- Meets IP54 Waterproof and dustproof rating, connectors of the meter protected by rubber cap, ensuring quality, reliability and service life especially in tough environments.

1.3 Features in pH Measurement (Applicable Models: PH800, PC800)

- 1 to 3 points auto calibration with calibration guide and auto-check function .
- Automatic recognition of up to 8 types of calibration solutions. Two series of standard solutions for choice: USA and NIST.
- Automatic display of electrode slope.

1.4 Features in Conductivity Measurement (Applicable Models: EC800, PC800)

- 1 to 4 points auto calibration with calibration guide and auto-check function.
- Single-tap switch among conductivity, TDS, salinity, and resistivity.

2. What's in the Box?

	Content	Quantity	PH800	EC800	PC800
1	PH800 Lab pH Meter	1	√		
2	EC800 Lab Conductivity Meter	1		√	
3	PC800 Lab pH/Conductivity Meter	1			√
4	201T-F 3-in-1 Combination pH Electrode ATC	1	√		√
5	2301T-F Conductivity Electrode ATC (ATC, k=1.0)	1		√	√
6	pH Standard Buffer (4.00 pH, 7.00 pH, 10.01pH/50mL)	1 bottle each	√		√
7	Conductivity Standard Solution (84μS, 1413μS, 12.88mS/50mL)	1 bottle each		√	√
8	9V Power Adapter	1	√	√	√
9	Flexible Electrode Holder	1	√	√	√
10	PC-Link Software Disk	1	√	√	√
11	USB Cable	1	√	√	√
12	Instruction Manual	1	√	√	√
13	Quick Manual	1	√	√	√

3. Technical Specifications

pH	Measuring Range	(-2.00 ~ 19.99) pH	PH800 PC800
	Resolution	0.1/0.01 pH	
	Accuracy	±0.01 pH ±1 digit	
	Temp. Compensation Range	(0 ~ 100°C) Automatic or Manual	
mV	Measuring Range	±1999mV	
	Resolution	1mV	
	Accuracy	±0.1% FS ±1 digit	
Cond.	Measuring Range	Conductivity: (0~200)mS/cm, including 5 ranges: (0.00~19.99)µS/cm, (20.0~199.9)µS/cm, (200~1999)µS/cm, (2.00~19.99)mS/cm, (20.0~199.9)mS/cm.	EC800 PC800
		TDS: (0 ~ 100)g/L, including 5 ranges: (0.00~9.99)mg/L, (10.0~99.9)mg/L (100~999)mg/L, (1.00~9.99)g/L (10.0~99.9)g/L	
		Salinity: (0~100) ppt, including 2 ranges: (0 ~ 9.99) ppt, (10.0 ~99.9) ppt	
		Resistivity: (0~100) MΩ·cm, including 6 ranges: (0 .0~ 99.9)Ω·cm, (100~ 999)Ω·cm (1.00~ 9.99)KΩ·cm, (10.0~ 99.9)KΩ·cm (100~ 999)KΩ·cm, (1.0~ 99.9)MΩ·cm	
	Resolution	Conductivity: 0.01/0.1/1µS/cm, 0.01/0.1 mS/cm TDS: 0.01/0.1/1mg/L, 0.01/0.1g/L Salinity: 0.01/0.1/1 ppt, 0.01/0.1 ppt. Resistivity: 0.1/1Ω·cm, 0.01/0.1/1 KΩ·cm, 0.1 MΩ·cm	
	Accuracy	±1.0% FS ±1 digit	

Cond.	Temp. Compensation Range	(0 ~ 50℃) Automatic or Manual	EC800 PC800
	Cell constant	0.1 / 1 / 10 cm ⁻¹	
Temp.	Measuring Range	0 to 100℃	PH800 EC800 PC800
	Resolution	0.1 °C	
	Accuracy	±0.5 °C ±1 digit	
Other	Data Storage	PH800/EC800: 500 Groups ; PC800: 1000 Groups	
	Storage Content	Numberings, Date, Time, Measurements, Unit, Temperature	
	Data Output	USB	
	IP Rating	IP54 Waterproof and Dustproof	
	Calibration	pH: 1~3 points , Conductivity: 1~4 points	

4. Instrument Description

4.1 LCD Display

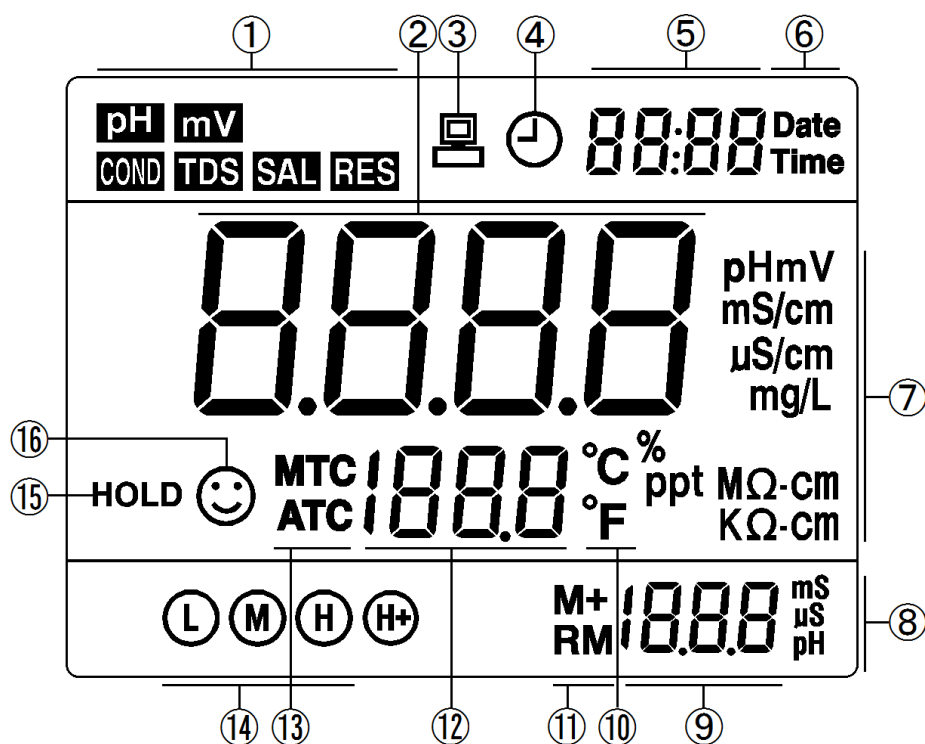


Diagram-1

- ① -- Measuring parameters
- ② -- Measuring value
- ③ -- USB data communication.

The meter is connected to your PC when this icon is displayed.

- ④ -- Timing storage
- ⑤ -- Date, time, and reminder icons
- ⑥ -- "Date" and "Time"
- ⑦ -- Measuring unit
- ⑧ -- Measuring unit in calibration
- ⑨ -- Calibration value, numberings of data storage, and reminder icons
- ⑩ -- Temperature unit
- ⑪ -- Data storage and display icon
M+– Data storage; RM– Recall saved data
- ⑫ -- Temperature and reminder icons
- ⑬ -- ATC–Auto Temperature Compensation;
MTC– Manual Temperature Compensation
- ⑭ -- Completed calibration icons
- ⑮ -- Auto-Lock reading mode
- ⑯ -- Stable reading icon

4.2 Keypad

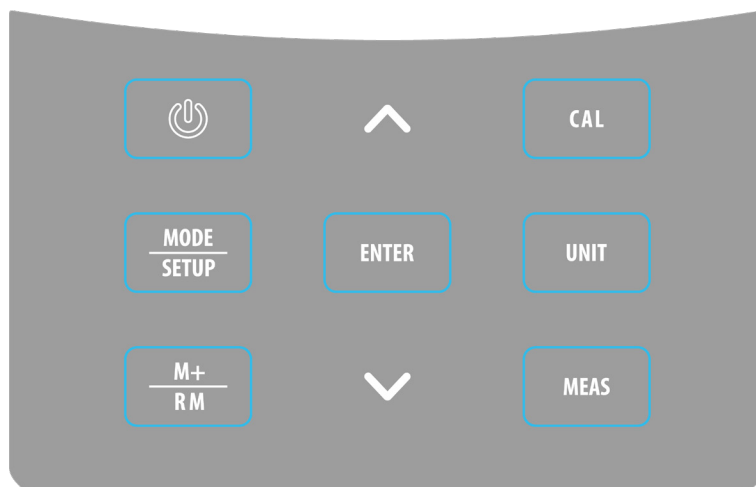


Diagram-2

Short Press – <1.5 s; Long Press – >1.5 s

Table -1 Keypad operations and functions

Keypad	Operation	Functions
	Short press	Power on/off
	Short press	Choose measuring mode: <ul style="list-style-type: none"> pH Meter: pH → mV pH/Conductivity Meter: pH → mV → COND
	Long press	Enter parameter setting
	Short Press	Enter calibration mode
	Short Press	<ul style="list-style-type: none"> In pH mode: choose resolution: 0.1 pH → 0.01pH In conductivity mode, choose measuring mode • COND → TDS → SAL → RES
	Short Press	<ul style="list-style-type: none"> In calibration mode: press to calibrate In parameter setting mode: press to confirm choice

	Short Press	Cancel any operation, the meter goes back to measurement mode
	Short Press	Save measuring data
	Long Press	Recall saved measuring data
 	Short Press or Long Press	<ul style="list-style-type: none"> • In manual temperature compensation (MTC) mode: Short press to adjust temperature, long press to adjust swiftly • In parameter setting mode: press to change the numbering of parameters in main menu and sub-menu. • In sub-menu, press to change parameters and settings. • In recall mode (RM), short press to change numberings, long press to change swiftly.

4.3 Connectors

Table-2 Connectors of Different Models

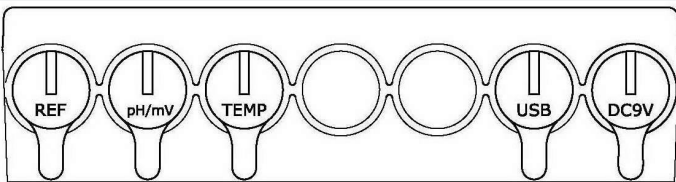
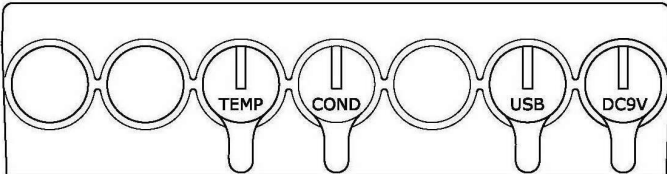
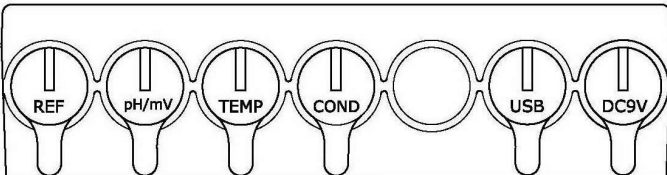
Connectors	Model
	PH800 pH Meter
	EC800 Conductivity Meter
	PC800 pH/Conductivity Meter

Table-3 Connector Name and Type

Icon	Connector Name	Connector Type
REF	Reference Electrode Plug	Φ2 Banana Plug
pH/mV	pH Electrode and ORP Electrode Plug	BNC Connector
TEMP	Temperature Plug (shared by pH and conductivity)	RCA Connector
COND	Conductivity Plug	BNC Connector
USB	USB Communication Plug	Standard Four Core Connector
DC9V	DC9V Power Plug	Φ2.5 Power Plug

special notes: don't pull out the power plug when the meter is working. It should be pulled out only after the meter is turned off.

4.4 Display Mode

4.4.1 Stable Reading Display Mode

When the measuring value is stable, the screen displays ☺ as shown in Diagram-3. If ☺ does not appear or is flashing, that means the measuring value has not been stable. Users should wait for the smiley face to appear and stay before recording the readings or conducting calibrations.

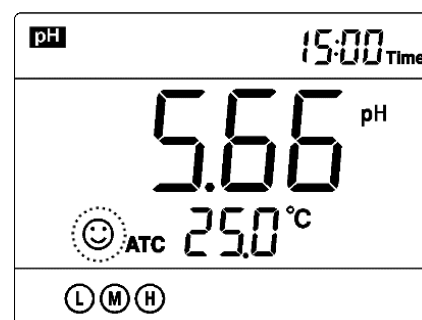


Diagram-3

4.4.2 Auto-Lock Display Mode

In parameter setting P5.4, select “On” to turn on the auto-lock mode, in which the measuring reading will be automatically locked after ☺ icon has been stable for 10 seconds, and the HOLD icon will come up as shown in Diagram-4. At this time, press **MEAS** to cancel the hold.

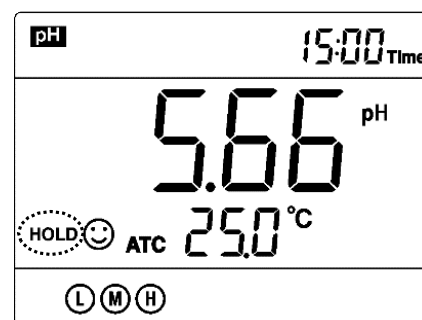


Diagram-4

4.5 Data Storage, Recall, and Deletion

4.5.1 Manual Data Storage

When reading is stable, short press , the screen will display M+ Icon along with the storage numbering, and the data will be saved, as shown in Diagram-5. For the storage capacity for each model, please refer to section 3.

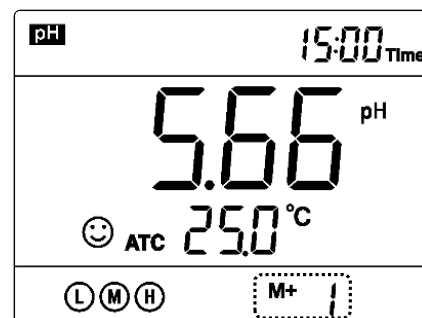


Diagram-5

4.5.2 Automatic Timing Data Storage

In parameter setting P5.1, set up the time for timing measurements

(e.g. 3 minutes).  be displayed, meaning the meter is ready to enter the auto timing storage mode. Short press , the  icon will start flashing, and the 1st group of data will be stored. 3 minutes later, the 2nd one will be stored. Diagram-6 shows that 8 groups of data have been automatically

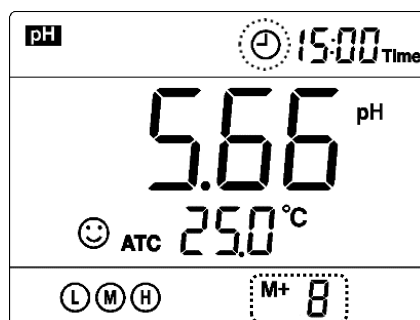


Diagram-6

stored. Short press  again, the  icon will stop flashing. The meter stops the auto timing storage. When in auto timing storage, manual storage is disabled. In parameter setting P5.1, set the time to 0 to exit auto timing storage.

4.5.3 Recall the Stored Data

In measurement mode, long press  (>1.5 s), the meter will recall the Measured value that was lastly saved, as shown in Diagram-7, displaying RM Icon and the numbering. Press  or  to recall other stored data. Long press  or  to change numberings swiftly.

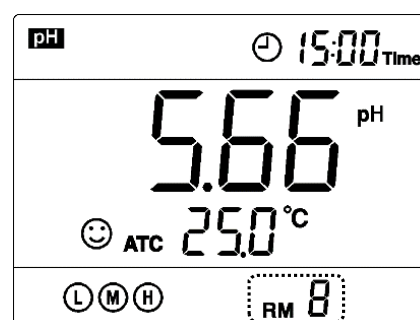


Diagram-7

4.5.4 Delete Stored Data

In parameter P5.3, select Yes to delete all the stored data. For details, please refer to section 7.5

5. pH Measurement

5.1 Preparation

5.1.1 Set up the Flexible Electrode Holder

The flexible electrode holder is composed of a base and an electrode holder. Place the electrode holder right above the metal stick on the base through the hole underneath; push down; and adjust the nut on electrode holder to finish installation.

5.1.2 Connect the power adapter

Plug the power adapter tightly into the DC9V socket. The power's voltage should meet the requirements listed on the power adapter.

5.1.3 Things needed in addition to what's included in the box

A clean cup, distilled or deionized water (8-16oz), and tissue papers for rinsing and drying the probe.

5.2 pH Electrode Information

5.2.1 pH Electrode

The meter comes with the 201T-F 3-in-1 combination electrode with a built-in temperature sensor, which enables the automatic temperature compensation. This electrode is only suitable for regular water solutions' pH testing. Please refer to Section 10 for ideal pH electrodes to use for other specific applications. The electrode's housing adopts polycarbonate materials, which is resistant to shocks and corrosion. A 3M KCL storage bottle comes with the electrode, which is for storing the electrode when not in use to keep the sensitivity of the probe.

5.2.2 Technical Specifications of the 201T-F pH Electrode

Measuring Range	0-14 pH
Temperature Range	32 to 176 °F (0 to 80 ° C)
Reference	Ag/AgCl
Junction	Ceramic
Connector	BNC/RCA
Temperature unit	30KΩ Thermistor
Electrode Dimension	Φ12×160 mm

5.2.3 Electrode Connectors

The pH electrode has two connectors: the BNC connector connects the pH probe; the RCA connector connects the temperature sensor. Plug these two connectors into ‘pH/mV” and “TEMP” sockets. Please note not to pull the cables in case of poor contact. Please keep the connectors clean and dry. Refer to section 4.7 regarding how to properly maintain the pH electrode.

5.2.4 How to Measure with the Electrode

Screw off the KCL storage bottle, and put it aside (do not dump or spill the KCL solution). Prepare a cup of distilled or deionized water, and rinse the electrode in it for a few seconds. Gently shake the probe to remove excess water, and dry it with clean tissue paper (DO NOT rub or wipe the probe, just use paper to dap off excess water). Gently stir the electrode for a few seconds after it’s dipped into the test solution and then let it stand still. Wait for the stable measurement (a smiley face appears and stays) and then record the readings. When test is finished, place the electrode in the storage bottle and tighten the cap in order to keep the sensitivity of the pH sensor.

Notes:

1) This default 201T-F electrode **will NOT** give accurate and stable pH readings when testing purified water such as distilled or deionized water. To measure distilled or deionized water's pH, users need to use a specialized pH electrode with the meter such as LabSen 803 Pure Water pH Electrode. Contact us at 1-614-285-3080 or info@aperainst.com for more details. When testing purified water like spring water or drinking water, it will take longer for the readings to get stabilized (typically 3-5 minutes) because there is very few ions left to be detected by the sensor in those purified water.

2) Purified water such as distilled water and deionized water are recommended **for only** rinsing the probe for the best result.

3) The 201T-F pH Electrode is **NOT** suitable for testing high-temperature samples (>176°F). Testing high-temperature samples could cause permanent damage to the electrode. A specialized electrode such as LabSen 213 (up to 225°F) should be used in this situation.

4) The electrode housing is polycarbonate. When using cleaning solutions, take cautions on carbon tetrachloride, trichlorethylene, tetrahydrofuran, acetone, etc., which will dissolve the housing and invalidate the electrode.

5.3 Information regarding pH Calibration

5.3.1 Standard Buffer Solutions

The instrument adopts two series of standard buffer solutions: USA and NIST as shown in Table-4. Users can select which one to use in P1.1 (refer to section 7.3)

Table-4 pH Standard Buffer Series

Calibration Icons		pH Standard Buffer Series	
		USA Series	NIST Series
3 Points Calibration		1.68 pH & 4.00 pH	1.68 pH & 4.01pH
		7.00 pH	6.86 pH
		10.01pH & 12.45 pH	9.18 pH & 12.45 pH

5.3.2 pH Calibration Modes

Users can adopt 1 to 3 points to calibrate as needed. In the 3-point calibration mode, the 1st point must be 7.00 pH (or 6.86 if using NIST). Then choose other calibration solutions to conduct 2nd and 3rd points (see Table-5 for details). In the process of calibration, the meter will display the electrode's slope in acid and alkaline ranges.

Table-5 3-point Calibration Mode

	USA series	NIST series	Calibration icon	When to adopt
1-point Calibration	7.00 pH	6.86 pH		Accuracy $\geq \pm 0.1$ pH
2-point Calibration	7.00 pH $\rightarrow$ 4.00 pH or 1.68 pH	6.86 pH $\rightarrow$ 4.01 pH or 1.68 pH		Measuring range: 0 to 7.00 pH
	7.00 pH $\rightarrow$ 10.01 pH or 12.45 pH	6.86 pH $\rightarrow$ 9.18 pH or 12.45 pH		Measuring range: >7.00 pH
3-point Calibration	7.00pH $\rightarrow$ 4.00 pH or 1.68 pH $\rightarrow$ 10.01pH or 12.45 pH	6.86pH $\rightarrow$ 4.01 pH or 1.68 pH $\rightarrow$ 9.18 pH or 12.45 pH		Wide measuring range

5.3.3 How often to calibrate

The frequency that you need to calibrate your meter depends on the tested samples, performance of electrodes, and the requirement of the accuracy. For High-Accuracy measurements ($\leq \pm 0.03\text{pH}$), the meter should be calibrated before test every time; For ordinary-accuracy Measurements ($\geq \pm 0.1\text{pH}$), once calibrated, the meter can be used for about a week or longer. In the following cases, the meter must be re-calibrated:

- The electrode hasn't been used for a long time or the electrode is brand new.
- After measuring strong acid ($\text{pH} < 2$) or strong base ($\text{pH} > 12$) solutions.
- After measuring fluoride-containing solution and strong organic solution
- There is a big difference between the temperature of the test sample and the temperature of the buffer solution that is used in the last calibration.

5.3.4 Calibration Reminder Function

Preset the interval between calibrations (starting from the time when you set it), and then the meter will remind you to calibrate at the end of that interval. For details, please see P1.2 (Section 7.3). When the (refer to section 7.3), preset time is reached, Er6 icon will be displayed at the lower right corner of the LCD (as showed in Diagram-8). At the time, the meter can still be operated. It is just reminding you to do calibration in order to ensure the accuracy. After calibration, the Er6 icon will disappear; To make it disappear, users can also choose "No" in P1.2 in parameter setting.

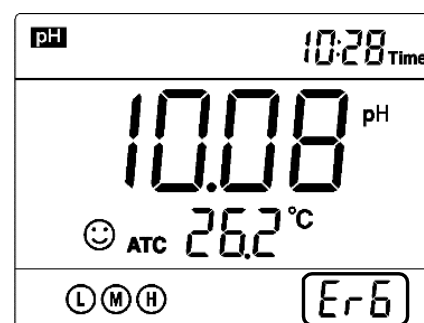


Diagram-8

5.3.5 Check calibration date

In this mode, users can see the date and time of last calibration to help determine if there is a need to re-calibrate. For details, please see parameter setting P1.3 (Section 7.3)


5.3.6 pH Isothermal Measurement Principle

According to the pH isothermal measurement principle, the closer the test sample's temperature is to the calibration solution's, the higher the accuracy of the measurement. This principle is recommended to follow when conducting tests for the best result. Example: If users were to test samples at 150°F, we recommend warming up the calibration solutions to the same temperature before performing calibrations in order to get the most accurate readings.

5.4 pH Calibration (Take 3-point calibration as an example)

5.4.1 Press **CAL** to enter calibration mode. CAL1 icon will flash in the upper right corner of the LCD. 7.00 pH will flash in the lower right corner of the LCD, reminding you to use pH 7.00 buffer solution to conduct 1st point of calibration.

5.4.2 Use distilled water to rinse off electrode and then dry it. Dip it into pH 7.00 buffer solution, stir gently and let it stand still and wait for the reading to become stable. In the lower right corner of LCD, the process of auto recognizing the buffer solution will be displayed. Pressing **ENTER** before the buffer is recognized will generate Er2 (please refer to table 6).

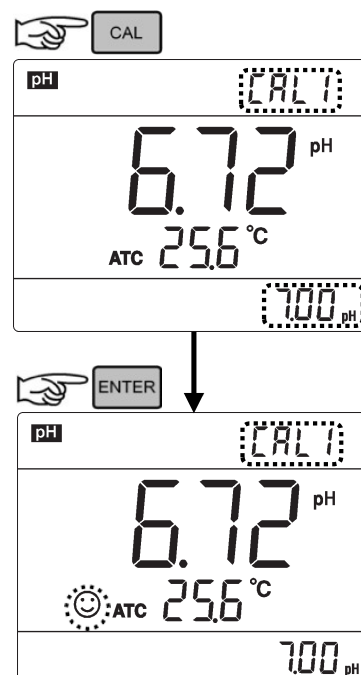

5.4.3 When the meter locks 7.00 pH, 😊 displays on LCD. Press **ENTER** key to calibrate the meter. **End** icon appears after calibration is done. The 1st point calibration is finished. In the meanwhile, CAL2 will flash at the upper right corner, and 4.00 pH & 10.01 pH will flash alternately at the bottom right, indicating using pH 4.00 or pH 10.01 buffer solution to make the 2nd point calibration.

5.4.4 Take out pH electrode, rinse it in distilled water, dry it, and dip it into pH 4.00 buffer solution. Stir the solution gently and let stand still in the buffer solution until a stable reading is reached. The meter's display will show the recognition process of calibration buffer solution at the lower right of LCD. When the meter recognizes 4.00 pH, 😊 displays on LCD. Press **ENTER** key to calibrate the meter. End icon and electrode slope of acidity range display after calibration is done. In the meanwhile, CAL3 will flash at the upper right corner of the LCD, and 10.01 pH will flash at the lower right, indicating using pH 10.01 buffer solution to make the 3rd point calibration.

5.4.5 Take out pH electrode, rinse it in distilled water, dry it, and dip it into pH 10.01 buffer solution. Stir the solution gently and let it stand still in the buffer solution until a stable reading is reached. The meter's display will show recognition process of calibration buffer solution at the bottom right of LCD. When the meter


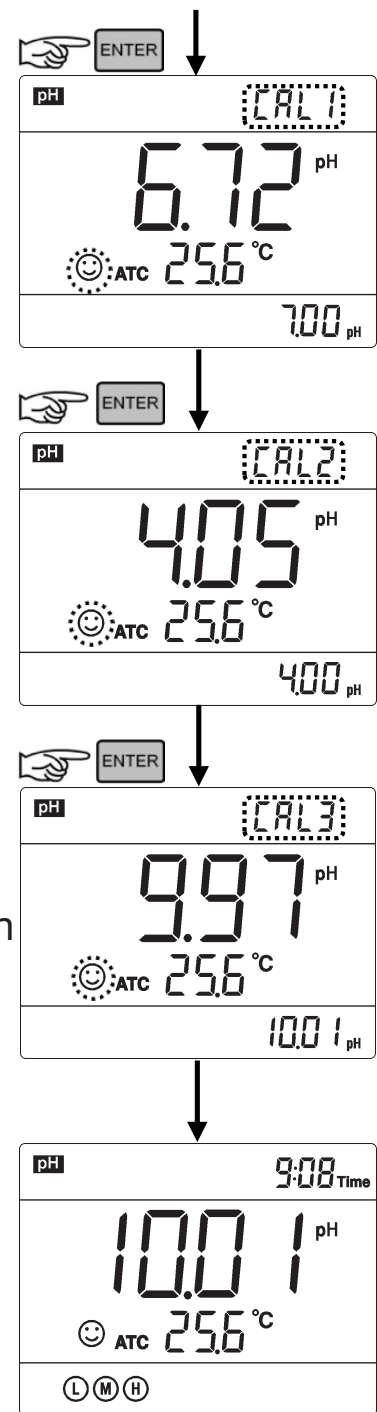


Diagram-9

recognizes 10.01 pH, 😊 displays on LCD. Press  key to calibrate the meter. End icon and electrode slope of alkalinity range display after calibration is done. The meter returns to the measurement mode, displays stable measuring value and calibration guide icons    . Please see Diagram-9 for the above calibration process.

5.5 Sample Measurement

5.5.1 Rinse the pH electrode in distilled water or purified water, dry it, and dip it into sample solution. Stir the solution gently and let it stand still in the sample solution until 😊 icon appears and stays on LCD, get the pH reading, which is pH value of sample solution, please refer to Diagram-10 for calibration and measurement process of the pH meter.

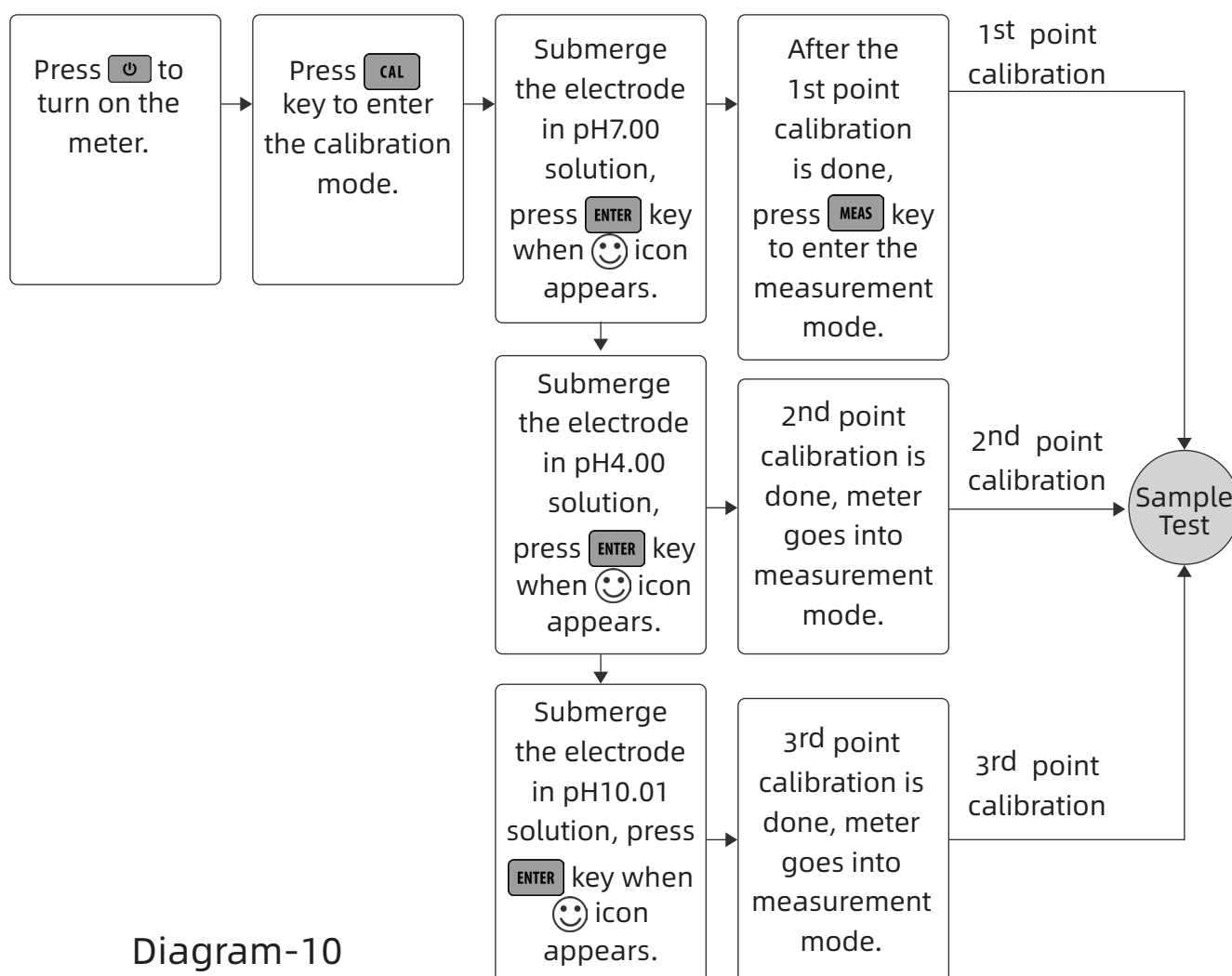

5.5.2 Self-Diagnosis Information

Among the calibration and measurement procedure, the meter has self-diagnosis function, reminding related information. Please refer to table-6.

Table – 6 Self-diagnosis information

Display Icons	Contents	Checking
<i>Er 1</i>	Wrong pH buffer solution or the buffer solution out of range.	<ol style="list-style-type: none"> 1. Check whether pH buffer solution is correct. 2. Check whether the meter connects the electrode properly. 3. Check whether the electrode is damaged.
<i>Er 2</i>	Press ENTER when reading is not stable during calibration.	Press ENTER key when 😊 icon appears and stays.
<i>Er 3</i>	During calibration, readings being unstable for over 3 minutes	<ol style="list-style-type: none"> 1. Check whether there are bubbles in glass bulb. 2. Replace with a new pH electrode.
<i>Er 4</i>	pH electrode zero electric potential out of range ($<-60\text{mV}$ or $>60\text{mV}$)	<ol style="list-style-type: none"> 1. Check whether there are bubbles in glass bulb. 2. Check whether pH buffer solution is correct.
<i>Er 5</i>	pH electrode slope out of range ($<85\%$ or $>110\%$)	<ol style="list-style-type: none"> 3. Check whether pH electrode is damaged. 4. The pH electrode is aged. Replace a new one
<i>Er 6</i>	The pre-set calibration due date is reached. Time to do perform a new calibration	Press CAL key to perform calibration or cancel due calibration setup in parameter P1.2.

* If you find any air bubble in the glass bulb of the pH sensor, simply shake the probe for a few times to remove it. The existence of an air bubble in the glass bulb will significantly decrease the stableness of measurement.

5.5.3 The instrument has a function to return to factory default setting, which can be set up in P1.4 (refer to section 7.3). Returning to factory default setting is to restore the meter to theoretical value (zero potential pH is 7.00, slope is 100%), and set all the parameters to default settings (see appendix 1). When the meter's calibration or measurement is performing abnormally, users can use this function to let the meter return to factory default mode, and then conduct calibration and test again. Please note that this function is irreversible once used.

5.6 Maintenance of the pH Electrode

5.6.1 Daily maintenance

The soaking solution contained in the supplied protective bottle is used to maintain activation in the glass bulb and junction. Loosen the capsule, remove the electrode and rinse in distilled water before taking a measurement. Insert the electrode and tighten the capsule after measurements to prevent the solution from leaking. If the soaking solution is turbid or moldy, replace the KCL solution (SKU: AI1107). Note: using other brands' storage solutions could cause permanent damage to the pH electrode because different chemicals are used. The electrode should not be soaked in distilled water, protein solution or acid fluoride solution for long periods of time. In addition, do not soak the electrode in organic lipids.

5.6.2 Calibration buffer solution

For calibration accuracy, the pH of the standard buffer solution must be reliable. The buffer solution should be refreshed often, especially after heavy use. We recommend 10-15 times of use before replacing the pH buffers.

5.6.3 Protect the glass bulb sensor

The sensitive glass bulb at the front of the combination electrode should not come in contact with hard surfaces. Scratches or cracks on the electrode will cause inaccurate readings. Before and after each measurement, the electrode should be washed with distilled water and dried. If a sample sticks to the electrode or it's contaminated, the electrode should be thoroughly cleaned using a soft brush and then rinsed with distilled water. After that, soak it in the KCL solution again for 6 hours.

5.6.4 Renew glass bulb

Electrodes that have been used over a long period of time will become aged. Submerge the electrode in 0.1mol/L hydrochloric acid for 24 hours, then wash the electrode in distilled water, then submerge it in KCL soaking solution for 24 hours. The method to prepare 0.1mol/L hydrochloric acid: dilute 9mL hydrochloric acid in distilled water to 1000mL. For serious passivation, submerge the bulb in 4% HF (hydrofluoric acid) for 3-5 seconds, and wash it in distilled water, then submerge it in the soaking solution for 24 hours to renew it.

5.6.5 Clean contaminated glass bulb and junction (please refer to Table-7)

Table-7 Clean contaminated glass bulb and junction

Contamination	Cleaning Solutions
Inorganic metal oxide	Dilute acid less than 1mol/L
Organic lipid	Dilute detergent (weak alkaline)
Resin macromolecule	Dilute alcohol, acetone, ether
Proteinic haematocyte sediment	Acidic enzymatic solution (saccharated yeast tablets)
Paints	Dilute bleach, peroxide

6. mV Measurement

6.1 mV Measurement

Press  to open, and press  to switch the meter to **mV** mode. mV mode is a mode for ORP and Ion electrode measurement. Connect ORP electrode and dip it in sample solution, stir the solution briefly and allow it to stay in the solution until  icon appears and get the reading. ORP means Oxidation Reduction Potential. The unit is mV.

6.2 Notes for ORP measurement

6.2.1 An ORP electrode is needed for ORP measurement (e.g. Apera 301Pt-C). ORP measurement does not require calibration. When the user is not sure about ORP electrode quality or measuring value, use ORP standard solution to test mV value and see whether ORP electrode or meter works properly. Table-8 is the data of standard ORP solution for 222 mV.

Table-8 ORP Standard Buffer Solution (222mV±15mV, 25°C)

°C	10	15	20	25	30	35	38	40
mV	242	235	227	222	215	209	205	201

6.2.2 Clean and activate ORP electrode

After the electrode has been used over a long period of time, the platinum surface will get polluted which causes inaccurate measurement and slow response. Please refer to the following methods to clean and activate ORP electrode:

- For inorganic pollutant, submerge the electrode in 0.1mol/L dilute hydrochloric acid for 30 minutes, then wash it in purified water, then submerge it in the soaking solution for 6 hours.
- For organic or lipid pollutant, clean the platinum surface with

detergent, then wash it in purified water, then submerge it in the soaking solution for 6 hours.

- c) For heavily polluted platinum surface on which there is oxidation film, polish the platinum surface with toothpaste, then wash it in purified water, then submerge it in the soaking solution for 6 hours.

7. Conductivity Measurement

7.1 Set up flexible Electrode Holder (refer to 4.1)

7.2 Information regarding the Conductivity Electrode

7.2.1 Conductivity electrode

EC800 and PC800 are equipped with model 2301T-F plastic conductivity electrode with constant $K=1.0$ and a built-in temperature sensor, which allows for automatic temperature compensation. When submerge the conductivity electrode in solution, stir the solution briefly to eliminate the air bubbles and improve response and stability.

7.2.2 Conductivity electrode constant

The meter matches conductivity electrodes of three constants: $K=0.1$, $K=1.0$ and $K=10.0$. Please refer to table-9 for measuring range. Set constant in parameter setting P2.1 and refer to clause 8.4

Table-9 Electrode constant and measuring range

Range	< 20 $\mu\text{S}/\text{cm}$	1.0 $\mu\text{S}/\text{cm}$ to 100 mS/cm			> 100 mS/cm
Conductivity electrode constant	K=0.1	K=1.0			K=10
Standard solution	84 $\mu\text{S}/\text{cm}$	84 $\mu\text{S}/\text{cm}$	1413 $\mu\text{S}/\text{cm}$	12.88 mS/cm	111.8 mS/cm
Electrode's model	DJS-0.1-F conductivity electrode	2301T-F conductivity electrode			2310T-F conductivity electrode

Note: When testing ultra-purified water with conductivity less than 1.0 $\mu\text{S}/\text{cm}$, a flow test should be conducted in a flow cell.

7.3 Information regarding Conductivity Calibration

7.3.1 Conductivity Standard Calibration Solutions

The meter uses conductivity standard solution of 84 $\mu\text{S}/\text{cm}$, 1413 $\mu\text{S}/\text{cm}$, 12.88 mS/cm and 111.8 mS/cm . The meter can recognize the standard solution automatically, and perform 1~4 points of calibration. The calibration indication icons correspond to the four standard values are in Table-10.

Table-10 Conductivity Calibration Solutions

Calibration Icons	Calibration Solutions
Ⓕ	84 $\mu\text{S}/\text{cm}$
Ⓜ	1413 $\mu\text{S}/\text{cm}$
ⓗ	12.88 mS/cm
ⓗ+	111.8 mS/cm

7.3.2 How often to calibrate

- a. The meter has been calibrated before leaving the factory and can generally be used right out of the box.
- b. Normally perform calibration once per month.
- c. For high accuracy measurements or large temperature deviation from the reference temperature (25° C), perform calibration once per week.
- d. Use conductivity standard solution to check whether there is error. Perform calibration if error is large.
- e. For new electrode or the meter has been set to factory default, perform 3-point or 4-point calibration. For general use, choose standard solutions that are closer to the sample solution to perform 1-point or 2-point calibration.

7.3.3 Single point and multi-point calibration

If 1-point calibration is conducted after 3-point or 4-point calibration being done, the previous calibration values in the same range will be replaced. In the meanwhile, the meter will display the 1-point calibration's icon, and the other calibration icons will be removed, but the chip will still store the data from the last calibration. When conducting multi-point calibration, users should follow the sequence from low conductivity to high in case the high concentration solution contaminating the low ones.

7.3.4 Reference Temperature

The factory default setting for reference temperature is 25 °C . The reference temperature can be set from 15°C to 30°C . Users can set it up in parameter setting P2.4 (see Section 7.4 for details).

7.3.5 Temperature compensation coefficient

The temperature compensation coefficient of the meter setting is 2.0%/°C . However, the conductivity temperature coefficient is different from solutions and concentration. Please refer to Table - 11 and the data collected during testing. Do the setting in P2.5. (see section 7.4 for more).

Note: When the coefficient for the temperature compensation is set to 0.00 (no compensation), the measurement value will be based on the current temperature.

Table-11 Temperature compensation coefficient of special solutions

Solution	Temperature compensation coefficient
NaCl solution	2.12 %/ °C
5% NaOH solution	1.72 %/ °C
Dilute ammonia solution	1.88 %/ °C
10% hydrochloric acid solution	1.32 %/ °C
5% sulfuric acid solution	0.96 %/ °C

7.3.6 Precaution for calibration solution's contamination

Conductivity standard solution has no buffer. Please avoid being contaminated during usage. Before submerging the electrode in standard solution, please wash the electrode and allow it to dry, especially for standard solution of low concentration (84μS/cm). The contaminated standard solution will affect accuracy.

7.3.7 Calibration Reminder Setup

Preset the interval between calibrations (starting from the time when you set it), and then the meter will remind you to calibrate at the end of that interval. For details, please see P2.3 (Section 7.4). When the preset time is reached, **Er6** icon will be displayed at the lower right corner of the LCD (as showed in Diagram-11). At the time, the meter can still be operated. It is just reminding you to do calibration in order to ensure the accuracy. After calibration, the **Er6** icon will disappear; To make it disappear, users can also choose “No” in P2.3 in parameter setting.

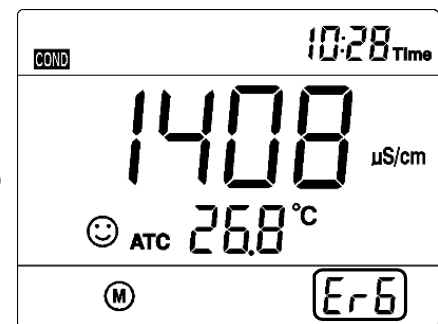


Diagram-11

7.3.8 Check calibration date

In this mode, users can see the date and time of last calibration to help determine if there is a need to re-calibrate. For details, please see parameter setting P2.4 (Section 7.4)

7.4 Conductivity Calibration (take 1413 $\mu\text{S}/\text{cm}$ as an example)

7.4.1 Rinse the electrode in purified water, allow it to dry, wash with a little of standard solution and submerge it in 1413 $\mu\text{S}/\text{cm}$ standard solution. Stir briefly and allow it to stay in the solution until a stable reading is reached.

7.4.2 Press **CAL** key to enter the calibration mode. The meter's display will show blinking—**CAL** at the top right, and scanning and locking process of calibration solution at the bottom right. **Er2** icon will appear if press **ENTER** before value is locked. See table-6

7.4.3 When the meter locks 1413 μS , stable icon ☺ will display on LCD. Press **ENTER** key to complete calibration. End will appear and the meter will return to measuring mode and **(M)** is displayed on bottom left of the LCD screen. See Diagram-12.

7.4.4 Press **MEAS** key before confirmation to exit calibration mode (calibration will not be completed).

7.4.5 For multi-point calibrations, please repeat the steps in 7.4.1 to 7.4.3 until all calibrations are finished. The meter can be calibrated in the same calibration solution repeatedly until the reading is stable and repeatable.

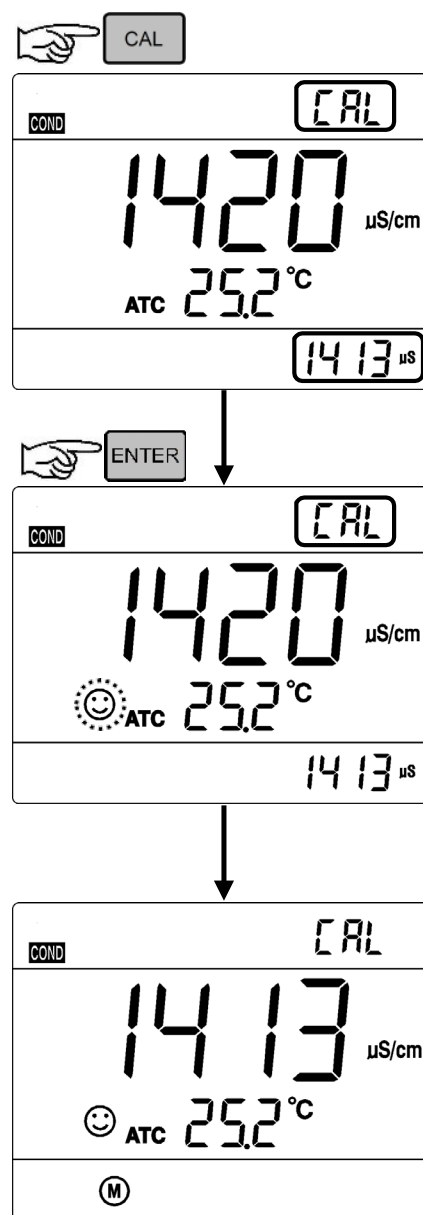


Diagram-12

7.5 TDS, Salinity, Resistivity & Conductivity

7.5.1 TDS and conductivity is linear related. The conversion factor is 0.40-1.00. Adjust the factor from parameter P2.7. The factory default setting is 0.71 and please refer to section 7.4. Salinity and resistivity are all proportional to conductivity. The calculating formula is preset in the meter. So the meter only needs to be calibrated in Conductivity mode, then after calibration of conductivity, the meter can switch from conductivity to TDS, salinity, or resistivity.

7.5.2 Adjust TDS conversion factor in parameter setting P2.7 according to the data collected during testing and experience.

Table - 12 lists some commonly used Conductivity and TDS conversion factors. This is for your reference only.

Table-12 Conductivity and TDS conversion factors

Conductivity of solution	TDS conversion factor
0~100 $\mu\text{S}/\text{cm}$	0.60
100~1000 $\mu\text{S}/\text{cm}$	0.71
1~10 mS/cm	0.81
10~100 mS/cm	0.94

***1000 ppm = 1 ppt**

7.6 Sample test

7.6.1 Rinse conductivity electrode in distilled or purified water, dry it, and submerge it in the sample solution. Stir the solution briefly and allow it to stay in the sample solution until a stable reading is reached and 😊 icon appears on LCD, then get the reading value, which is the conductivity value of the solution.

7.6.2 Press  key to switch to TDS, Salinity, and Resistivity.

7.6.3 During the process of calibration and measurement, the meter has self-diagnosis functions, indicating the relative information as below:

Table-13 Self-diagnosis information of conductivity measurement mode

Display Icons	Contents	Checking
<i>Er 1</i>	Wrong pH buffer solution or the buffer solution out of range.	<ol style="list-style-type: none"> 1. Check whether pH buffer solution is correct. 2. Check whether the meter connects the electrode properly. 3. Check whether the electrode is damaged.
<i>Er 2</i>	Press ENTER when measuring value is not stable during calibration.	Press ENTER key after 😊 icon appears and stays.
<i>Er 3</i>	During calibration, the measuring value being unstable for over 3 minutes	<ol style="list-style-type: none"> 1.Shake the electrode to eliminate bubbles in electrode head. 2.Replace with new conductivity electrode.
<i>Er 6</i>	Enter in pre-set due calibration to remind calibration	Press CAL key to perform calibration or cancel due calibration setup in parameter P2.2.

7.6.4 Factory default setting

For factory default setting, please refer to parameter setting P2.7 (Section 7.4). With this function, all calibration data is deleted and the meter restores to the theory value. Some functions restore to the original value (refer to appendix -1). When calibration or measurement fails, please restore the meter to factory default setting and then perform re-calibration or measurement. Please

note once set the factory default, all the data deleted will be irretrievable.

7.7 Maintenance of the Conductivity Electrode

7.7.1 Always keep the conductivity electrode clean. Before taking a measurement, rinse the electrode in distilled water. It is recommended to rinse it again in the sample solution. When submerge the electrode in solution, stir the solution briefly to eliminate air bubbles and allow it to make the measurement fast and stable.

7.7.2 At the top of the electrode there is a water sponge (see Diagram-13), when water sponge is dry up, please refill with distilled water to keep it wet (not too much). Keep in wet condition helps the electrode's fast response time. Electrode may react slowly if kept dry for a long period. Submerge it in 12.88Ms buffer solution for 5 to 10 minutes or in tap water for 1 to 2 hours to renew it. Use distilled water to clean electrode after use.

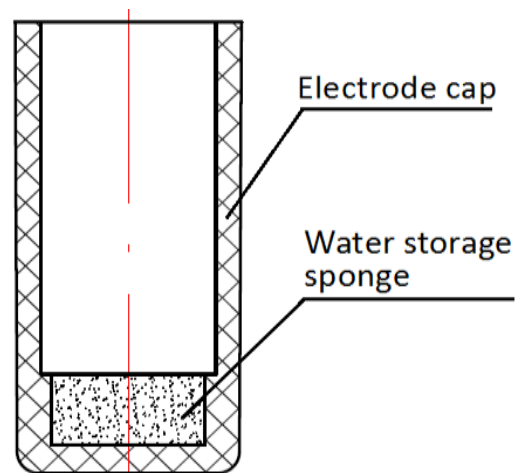


Diagram-13

7.7.3 The sensing rod of Model 2301T-F conductivity electrode is coated with platinum black to minimize electrode polarization and expand measuring range. The platinum black coating adopted our special processing technology, which improves the electrode performance and the firmness of the coating. If the platinum black electrode is stained, gently clean the electrode with soft brush in warm water containing detergent or alcohol.

8. Parameter Setting

8.1 Main menu

In the measurement mode, press  key to enter in P1.0. The meter has three main menus. Press  or  to switch main menu: P1.0 → P2.0 → P5.0. Please refer to Diagram - 14.

P1.0: pH parameter setting menu,
P2.0: Conductivity parameter setting menu,
P5.0: Basic parameter setting menu.

8.2 Sub-Menu

8.2.1 In P1.0, press  to enter the submenu P1.1 for pH setting, press  or  to change submenu: P1.1 → P1.2 → ... → P1.4. See 8.3 for details.

8.2.2 In P2.0, press  to enter the submenu P2.1 for conductivity setting, press  or  to change submenu: P2.1 → P2.2 → ... → P2.7. See 8.4 for details.

8.2.3 In P5.0, press  to enter the submenu P5.1 for basic parameter setting, press  or  to change submenu: P5.1 → P5.2 → ... → P5.6. See 8.5 for details.

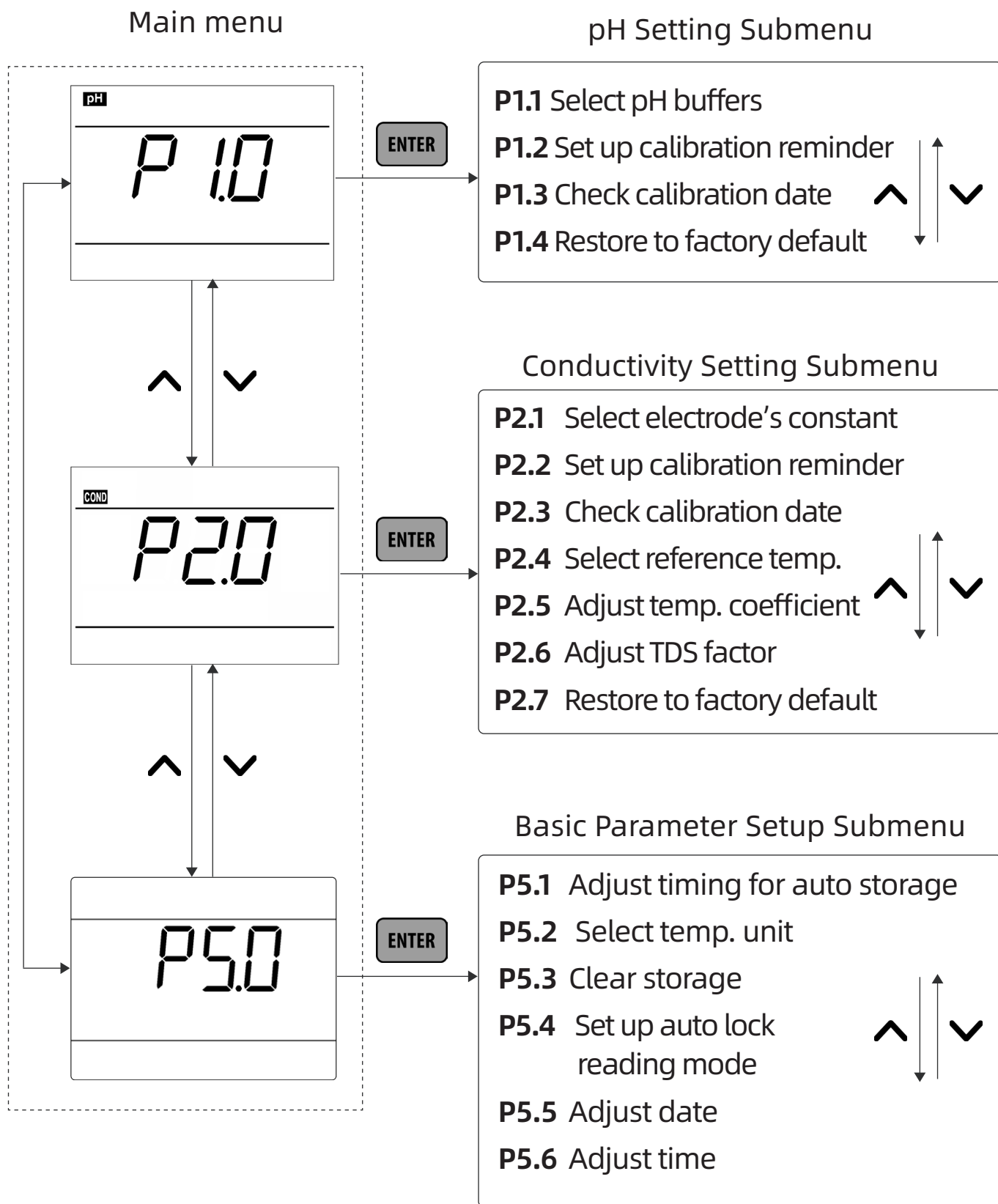
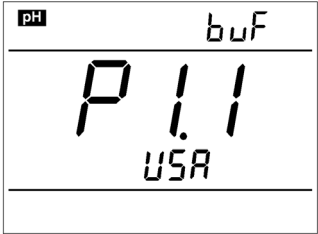
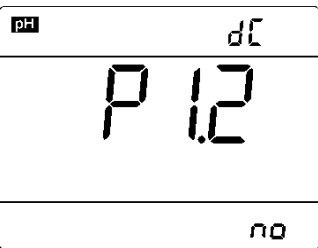
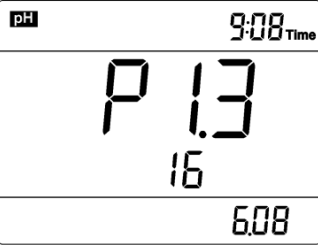
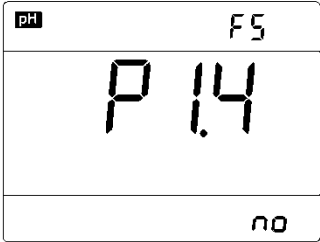




Diagram-14






8.3 pH Setting Sub-Menu (press or to switch)

	<p>P1.1 - Select pH buffer series (USA-NIST)</p> <ol style="list-style-type: none"> 1. In measurement mode, long press  to enter P1.0, press  to enter P1.1 2. Press , USA flashes, press  to choose USA → NIST, press  to confirm. USA–USA; NIS–NIST 3. Press  to enter P1.2, or press  to return to measurement mode.
	<p>P1.2 - Set up calibration reminder (No–H00–D00)</p> <ol style="list-style-type: none"> 1. Press , "No" flashes, press  to choose No → H00 → D00; NO–no setting; H00–set hours(0~99); D00–set days(0~99). 2. When "H" flashes, press , "00" flashes, press  to adjust hours, press  to confirm; When "D" flashes, press , "00" flashes, press  to adjust days, press  to confirm; When "No" flashes, press  to confirm. 3. Press  to enter P1.3, or press  to return to measurement mode.
	<p>P1.3 - Check Calibration Time</p> <ol style="list-style-type: none"> 1. The data and time of last calibration displayed on the left image is: On June 8, 2016 to 9:08 2. Press  to enter P1.4, or press  to return to measurement mode.





	<p>P1.4 – Set the Factory Default (No–Yes)</p> <p>Press ENTER , “NO” flashes, press ▲ to select NO → Yes, press ENTER to confirm, instrument back to measurement mode.</p> <p>No – Not restore; Yes – Restore factory default</p>
---	---



8.4 Conductivity Setting Sub-Menu (press **▲** or **▼** to switch)

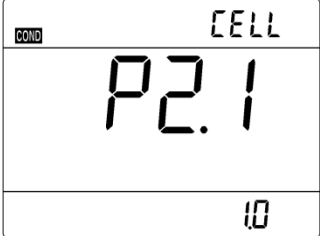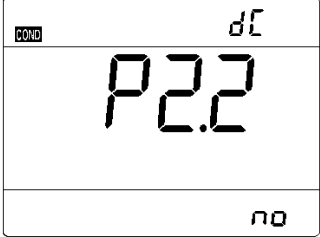	<p>P2.1 – Select electrode's constant (1.0–10.0–0.1)</p> <ol style="list-style-type: none"> 1. In P2.0, press to enter P2.1 2. Press ENTER , “1.0” flashes, press ▲ to choose 1.0 → 10.0 → 0.1; press ENTER to confirm. 3. Press ▲ to enter P2.2, or press MEAS key to return to measurement mode.
	<p>P2.2 – Set up calibration reminder (No–H00–D00)</p> <ol style="list-style-type: none"> 1. Press ENTER , “No” flashes, press ▲ to choose No → H00 → D00; NO–no setting; H00–set hours(0~99); D00–set days(0~99). 2. When “H” flashes, press ENTER , “00” flashes, press ▲ to adjust hours, press ENTER to confirm; When “D” flashes, press ENTER , “00” flashes, press ▲ to adjust days, press ENTER to confirm; When “No” flashes, press ENTER to confirm. 3. Press ▲ to enter P2.3, or press MEAS key to return to measurement mode.

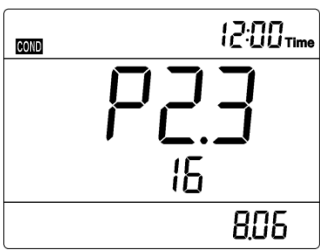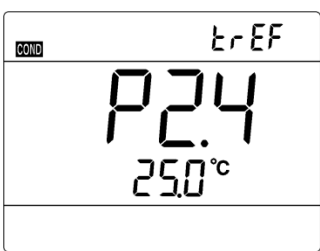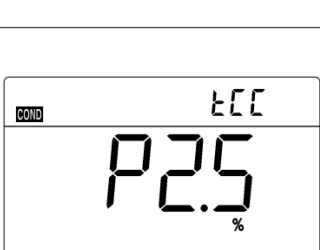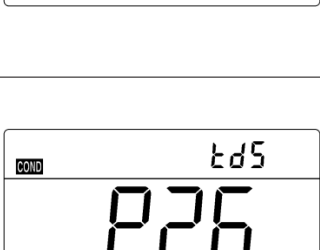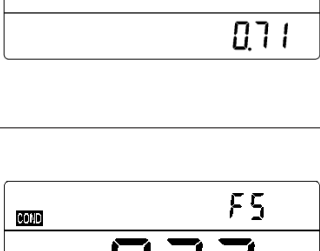	<p>P2.3 – Check Calibration Time</p> <ol style="list-style-type: none"> 1. The data and time of last calibration displayed on the left image is: On August 6, 2016 to 12:00. 2. Press  to enter P2.4, or press  key to return to measurement mode.
	<p>P2.4 – Select reference temperature (15.0°C–30.0°C)</p> <ol style="list-style-type: none"> 1. Press , “25.0°C” flashes, press  or  to adjust temp. from 15°C to 30°C ; press  to confirm. 2. Press  to enter P2.5, or press  key to return to measurement mode.
	<p>P2.5 – Adjust temp. compensation coefficient (0.00-9.99%)</p> <ol style="list-style-type: none"> 1. Press , “2.00” flashes, press  or  to adjust from 0.00-9.99; press  to confirm. 2. Press  to enter P2.6, or press  key to return to measurement mode.
	<p>P2.6 – Adjust TDS Factor (0.40-1.00)</p> <ol style="list-style-type: none"> 1. Press , “0.71” flashes, press  or  to adjust from 0.40-1.00; press  to confirm. 2. Press  to enter P2.7, or press  key to back to measurement mode.
	<p>P2.7 – Set the Factory Default (No–Yes)</p> <p>Press , “NO” flashes, press  to select NO → Yes, press  to confirm, instrument back to measurement mode.</p> <p>No – Not restore; Yes – Restore factory default</p>

8.5 Basic Parameter Setting Sub-Menu (press or to switch)

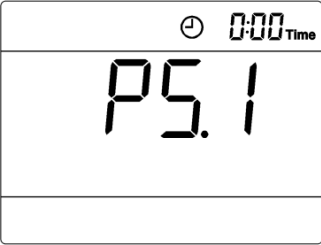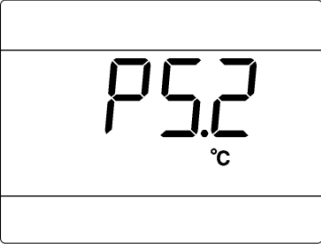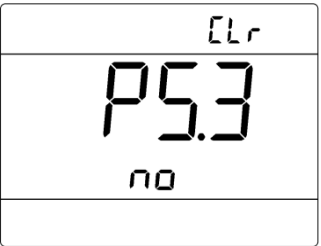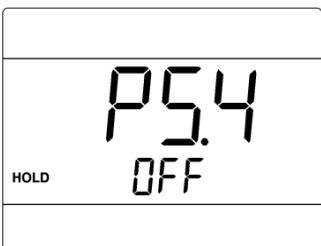	<p>P5.1 – Adjust timing for auto storage</p> <ol style="list-style-type: none"> 1. In P5.0, press  to enter P5.1 as shown in the left Diagram. 2. Press , ":00" flashes, press  or  to set minute (0-59), press  again, "0:" flashes, press  or  to set hours (0-99), press  to confirm. 3. Press  to enter P5.2, or press  key to return to measurement mode.
	<p>P5.2 – Choose Temp. Parameter (° C–° F)</p> <ol style="list-style-type: none"> 1. Press , "° C" flashes, press  to choose ° C → ° F, press  to confirm. 2. Press  to enter P5.3, or press  key to return to measurement mode.
	<p>P5.3 – Clear storage (No–Yes)</p> <ol style="list-style-type: none"> 1. Press , "No" flashes, press  to choose No → Yes; press  to confirm. No– not to clear data storage; Yes–clear all data storage 2. Press  to enter P5.4, or press  key to return to measurement mode.
	<p>P5.4 – Set up Auto-Lock reading mode (Off–On)</p> <ol style="list-style-type: none"> 1. Press , "Of" flashes, press  to choose Off → On, press  to confirm. Off– not set up, On–set up 2. Press  to enter P5.5, or press  key to return to measurement mode.

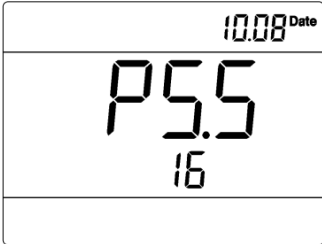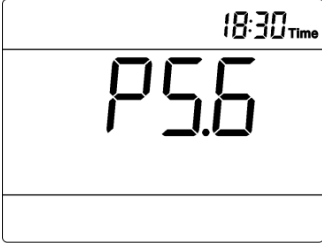	<p>P5.5 – Adjust Date</p> <ol style="list-style-type: none"> 1. Press ENTER to choose “Month”→ “Day”→ “Year”, press ▲ or ▼ to adjust date, press ENTER to confirm. 2. Press ▲ to enter P5.6, or press MEAS key to return to measurement mode.
	<p>P5.6 – Adjust Time</p> <ol style="list-style-type: none"> 1. Press ENTER to choose – “hour” → “minute”, press ▲ or ▼ to adjust time, press ENTER to confirm. 2. Press MEAS key to return to measurement mode.

9. USB Data Communication

9.1 Software Requirement

Two versions of PC-Link software are available: Web PC-Link-480 Software and the traditional PC-Link software.

Web PC-Link-480 Software is installation-free and supports both Windows and Mac computers. To use the Web PC-Link Software, open <https://pclink480.aperainst.io/PClink%20480%20V5> using Chrome, Edge, or Firefox, then connect the meter to the computer using a USB communication cable to start using the software.



Scan the code to access the
Web PC-Link webpage

If using the traditional PC-link software, the recommended requirement for the computer is a Windows-based Personal computer (with Microsoft Excel) with USB communication port. The following content is about the traditional PC-link software.

9.2 Software Interface

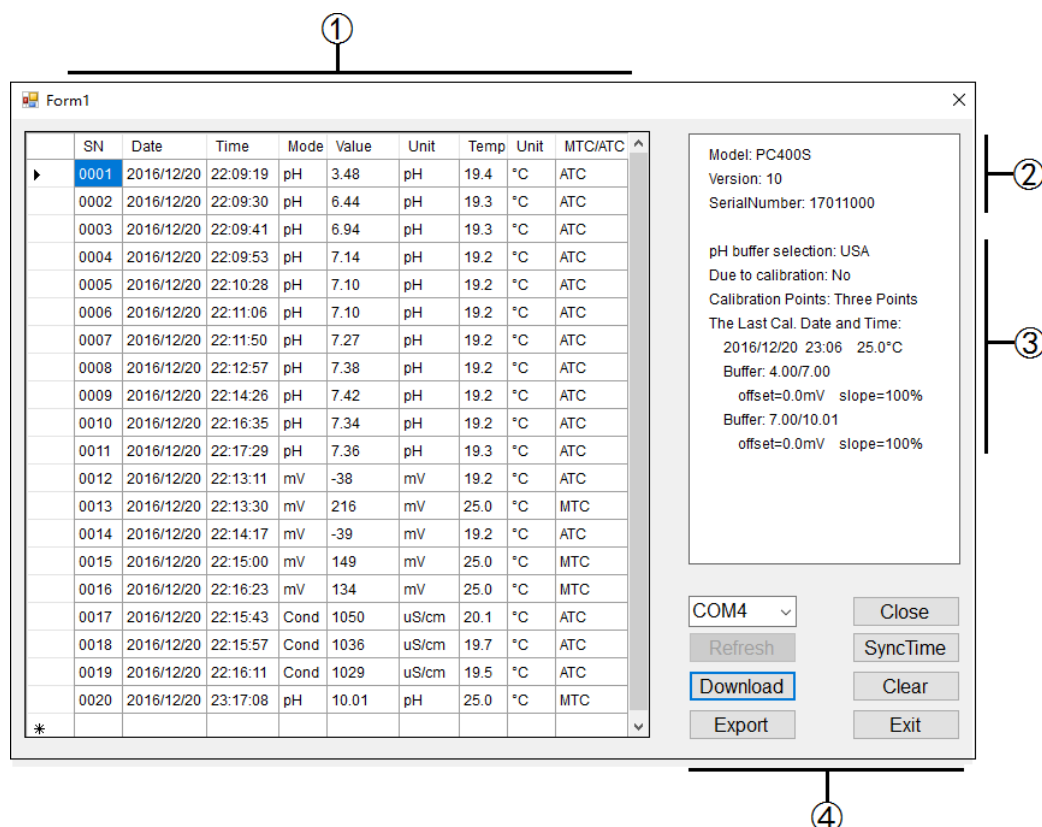


Diagram-15

① – Stored data

- Press “Download” Key to upload the data in the meter to the software, including date, time, measurements, temperature, and temperature compensation mode. The program will categorize the data by pH, mV, and Cond.
- Press  after the meter is connected to the computer, or set up auto timing storage function. In this case, all the measuring data will be uploaded to the software, and will not be stored in the meter. See Diagram-15.

- ② – Model and Numberings
- ③ – Calibration information
- ④ – COM Port #and other buttons

COM Port # is the number for the software to connect with the computer. Every computer's COM port is different. In Diagram-15, the COM port # is COM4.

- **Open/Close** – Click it to turn on/off the program. When turned on, the LCD will display 
- **Refresh** – COM port reset button; Click it to reset the COM port to COM1.
- **Sync Time** – Click it to sync the computer's time and date to the meter's.
- **Download** – Click it to upload the data stored in the meter to the computer.
- **Clear** – Click it to clear all the data.
- **Export** – Click it to export all the data to a Microsoft Excel file.
- **Exit** – Click it to exit the program.

9.3 Install the software

The PC-Link software works for all Windows based system (does not work for Mac). Insert the PC-Link USB flash drive into the computer, open the PC-Link folder where you will find the folder for PC-Link software and a zipped file for drivers. Typically, users can directly open the file with PCLink icon to use the software. If the meter cannot be connected to the computer, please install the USB driver (in the zipped file) before using the software.

9.4 Choose COM Port

Connect the meter to your computer with the USB cable. Open the

PC-Link program. click the arrow icon next to the COM1 Port, and click the bottom port number, and then click Open. The LCD will display  icon. If port number is hard to be confirmed, users can identify it in Windows' Device Manager.

9.5 Run Software

9.5.1 Upload stored data

Click Download key to upload the data stored in the meter to the software, including date, time, measurements, temperature, and temperature compensation mode. The program will categorize the data by pH, mV, and Cond.

9.5.2 Real-time storage

(a) When the program is running and meter is connected to computer, press  on the meter or set up auto timing storage to upload all the measuring data to the software. The data will not be stored in the meter in this mode.

(b) The mode and unit for real-time storage is the same as it is in the meter. Press  or  to make changes.

9.5.3 Data Processing

Click "Export" to export all the data to a Microsoft Excel file. Users can process, analyze, and print the data in the Excel file.

10. Recommended pH Electrodes for Different Applications

Application	Ideal Apera pH Electrodes to Use
General water solutions	LabSen 213
Dirty liquids such as wastewater, emulsions, suspensions, slurries, etc.	LabSen 333
Biotechnology and pharmaceuticals	LabSen 823
High-Temperature solutions	LabSen 863
Low-temperature liquid	LabSen 881
High salinity/ Strong alkaline solutions	LabSen 843
Complex and caustic solutions (e.g. electro-plating solutions)	LabSen 863
Pure Water (distilled/deionized/RO/ drinking/well/storm/boiler water)	LabSen 803
Liquid food and beverage (milk, yogurt, sauce, jam, etc.)	LabSen 823
Micro-volume liquid samples	LabSen 241-3, LabSen 241-3SP
Semi-micro volume liquid samples	LabSen 243-6
Frozen/raw meats	LabSen 763
Soil direct test	LabSen 553
Solid or semi-solid samples (cheese, dough, fruits, sushi rice, meat products, etc.)	LabSen 753
Strong acid solutions, HF containing solutions (HF concentration < 2000 ppm)	LabSen 833
Surface test (skin, paper, culture medium...)	LabSen 373
TRIS buffer solutions	LabSen 213/823

Application	Ideal Apera pH Electrodes to Use
Viscous liquid e.g. skincare products, coatings, syrups, sugarcane juice, etc.	LabSen 853-S/853-H
Organic solvents, non-aqueous solutions	LabSen 873

Note:

- Every pH electrode gradually ages and will eventually fail. A typical service life of a pH electrode is 1-2 years depending on many factors such as frequency of use, nature of test samples, how well it is maintained, etc. We recommend replacing your pH electrode at least every 1 to 2 years to guarantee the optimal performance.
- For more information, go to <https://aperainst.com/sensors/labсен> or contact us at info@aperaint.com or +1 (614) 285-3080.

11. Recommended Conductivity Electrodes for Different Applications

Application	Ideal Apera pH Electrodes to Use
General water solutions (K=1, 0 to 200 mS/cm)	2301T-F
High accuracy test for general water solutions (K=1, 0 to 200 mS/cm)	2401T-F
Pure water (K=0.1, 0 to 200 μ S/cm)	DJS-0.1-F

- For more information, go to <https://aperainst.com/sensors> or contact us at info@aperaint.com or +1 (614) 285-3080.

12. Replacement Parts and Accessories

WARNING



Use of non-approved parts may cause damage to the instrument or equipment malfunction. The replacement parts in this section are approved by the manufacturer.

Replacement Parts

Description	SKU	Quantity
LabSen211 glass pH combination electrode	AI3101	1 pc
2301T-F Conductivity Electrode ATC(ATC, k=1.0)	AI1305	1 pc
8 oz. pH/Conductivity Calibration Solution Kit (pH 7.00, 4.00, 10.01, 1413 μ S/cm, 12.88 mS/cm, 84 μ S/cm)+ 4oz. 3 M KCL storage solution + CalPod	AI1112	1 pc
602 flexible electrode holder	AI2810	1 pc
9V Power Adapter (120V Input) with US Plug	AI7108	1 pc
USB Cable for Data Transfer and Charging	AI2517	1 pc
Electrode Protective Cap	AI9201	1 pc
Small Electrode Storage Vial	AI2511	3 pcs

Accessories

Description	SKU	Quantity
301Pt-C ORP Electrode	AI1303	1 pc
222mV ORP Calibration Solution	AI1160	8 oz
111.8 mS/cm Conductivity Standard Calibration Solution	AI1134	4/8/16 oz
Electrode Cleaning Solution	AI1166	8 oz.
pH Electrode Maintenance Kit	AI1170	1 pc
9V Switching Power Adapter (Input: 100-240V) US Plug	AI7109	1 pc

Description	SKU	Quantity
9V Switching Power Adapter (Input: 100-240V) with Universal Plugs	AI7110	1 pc

For pH/Conductivity electrodes for special applications, refer to Section 10 & 11

13. Appendix I : Table of Parameter setting and factory default setting

Modes	Prompts	Parameter setting items	Abbreviation	Description	Restore to factory default setting
P1.0 pH	P1.1	Select Buffer Series	buF	USA-NIST	USA
	P1.2	Set calibration reminder	dC	No - H00 - D00	No
	P1.3	Check calibration date and time	/	-	-
	P1.4	Restore factory default	CAL	No - Yes	No
P2.0 Cond.	P2.1	Select electrode's constant	CELL	1.0 - 10.0 - 0.01	1.0
	P2.2	Set calibration reminder	dC	No - H00 - D00	No
	P2.3	Check calibration date and time	/	-	-
	P2.4	Select reference temperature	tREF	15~30°C	25°C
	P2.5	Adjust temperature compensation coefficient	tCC	0.00~9.99	2.00%
	P2.6	Adjust TDS factor	tDS	0.40~1.00	0.71
	P2.7	Restore factory default	FS	No - Yes	No

Modes	Prompts	Parameter setting items	Abbreviation	Description	Restore to factory default setting
P5.0 Basic Parameter	P5.1	Adjust timing for auto storage	/	-	-
	P5.2	Select temperature unit	/	°C - °F	-
	P5.3	Clear data storage	CLr	No - Yes	-
	P5.4	Set up auto-lock reading mode	/	Off-On	-
	P5.5	Adjust date	/	-	-
	P5.6	Adjust time	/	-	-

14. Appendix II: Icons and Abbreviation

Mode	Symbol	Abbreviation	Stand for
P1.0 pH	P1.1	buF	Standard buffers
	P1.2	dC	Due Calibration
	P1.3	/	/
	P1.4	FS	Factory default setting



Mode	Symbol	Abbreviation	Stand for
P2.0 Conductivity	P2.1	CELL	Cell
	P2.2	dC	Due Calibration
	P2.3	/	
	P2.4	tREF	Reference temperature
	P2.5	tCC	Temperature compensation coefficient
	P2.6	tDS	TDS coefficient
	P2.7	FS	Factory default setting
P3.0 Basic Parameters	P5.1	/	
	P5.2	/	
	P5.3	CLr	Clear readings
	P5.4	/	
	P5.5	/	
	P5.6	/	
Other		USA	United States of America
		n 15	Nist
		OFF	Off
		On	On
		no	No
		YES	Yes

15. Appendix III : Table of Self-diagnosis Symbol

Symbol	Self-Diagnosis Information	pH	Conductivity
<i>Er 1</i>	Wrong calibration solution or the meter recognition of calibration solution out of range.	√	√
<i>Er 2</i>	Press ENTER key when measuring value is not stable during calibration.	√	√
<i>Er 3</i>	During calibration, the measuring value is not stable for ≥ 3 min.	√	√
<i>Er 4</i>	pH electrode zero electric potential out of range (<-60mV or >60mV)	√	
<i>Er 5</i>	pH electrode slope out of range (<85% or >110%)	√	
<i>Er 6</i>	Enter in pre-set due calibration to remind re-calibration	√	√

16. Warranty

We warrant this instrument to be free from defects in material and workmanship and agree to repair or replace free of charge, at the option of APERA INSTRUMENTS, Co., Ltd, any malfunctioned or damaged product attributable to the responsibility of APERA INSTRUMENTS, Co., Ltd for a period of **THREE YEARS for the instrument and SIX MONTHS for the electrodes from the delivery.**

This limited warranty does NOT cover any issues due to:

- accidental damage
- transportation
- storage
- improper use
- failure to follow the product instructions or to perform any preventive maintenance
- materials, processes, systems or other matter not provided or authorized in writing by Apera
- unauthorized repair
- normal wear and tear
- external causes such as accidents, abuse, or other actions or events beyond our reasonable control

APERAINSTRUMENTS, LLC

Address: 6540 Singletree Dr., Columbus, Ohio 43229

Tel: 1-614-285-3080

Email: info@aperainst.com

Website: aperainst.com